

Anti-Mouse CD48-154Sm

Catalog: 3154004B

Package Size: 100 tests

Storage: Store product at 4°C. Do not freeze.

Reactivity: Mouse,

Clone: HM48-1


Isotype: Hamster IgG

Formulation: Antibody stabilizer with 0.05% Sodium Azide

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF® analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 µl for up to 3 X 10⁶ live cells in 100 µl. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

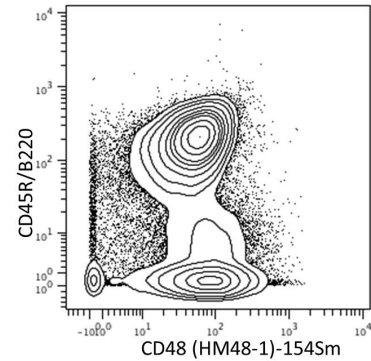

C57BL/6 mouse splenocytes stained with 176Yb-anti-CD45R/B220(RA36B2) and 154Sm-anti-CD48 (HM48-1).

Description

CD48, also known as BCM1, is a 45 kDa member of the signaling lymphocyte activation molecule (SLAM) family of proteins. It is expressed on T cells, B cells and monocytes/macrophages. CD48 plays a critical role in adhesion and T cell activation. It can act as a ligand for both CD244 and CD2 and, depending on the context, either inhibits or stimulates IFN γ production from the target cell. CD48 can be used in combination with CD150 and CD41 to identify CD150+CD48-CD41- hematopoietic stem cells (HSC).

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. **Highly Multiparametric Analysis by Mass Cytometry.** *J Immunol Methods* 361 (1-2):1-20, 2010.

For technical support visit fluidigm.com/support

North America +1 650 266 6100 | Toll-free: +1 866 358 4354 in the US | support.northamerica@fluidigm.com **Europe** +33 1 60 92 42 40 | support.europe@fluidigm.com

China (excluding Hong Kong) +86 21 3255 8368 | techsupportchina@fluidigm.com **Japan** +81 3 3662 2150 | techsupportjapan@fluidigm.com

All other Asian countries +1 650 266 6100 | techsupportasia@fluidigm.com **Central and South America** +1 650 266 6100 | techsupportlatam@fluidigm.com

For Research Use Only. Not for use in diagnostic procedures.

This product contains antibodies manufactured by and sold under license from BioLegend and licensees thereof.

Information in this publication is subject to change without notice. **Safety data sheet information** fluidigm.com/sds **Patent and license information** fluidigm.com/legalnotices | Fluidigm, the Fluidigm logo, and CyTOF are trademarks or registered trademarks of Fluidigm Corporation in the United States and/or other countries. © 2015 Fluidigm Corporation. All rights reserved. 07/2015