

Anti- β -actin-154Sm

Pathologist-Verified Clone for Imaging Mass Cytometry™

Catalog: 3154021D

Package size and concentration: 25 μ g, 0.5 mg/mL

Storage: Store at 4 °C. Do not freeze.

Reactivity: Rat, Mouse, Human

Clone: 2F1-1

Isotype: Mouse IgG2b

Formulation: Antibody stabilizer with 0.05% sodium azide

Application: IMC-Paraffin

Technical Information


Application: The metal-tagged antibody is designed and formulated for the application of Imaging Mass Cytometry (IMC™) using the Fluidigm Hyperion™ Imaging System on formalin-fixed, paraffin-embedded (FFPE) tissue sections.

Quality control: Each lot of conjugated antibody is quality control-tested by Imaging Mass Cytometry on tissue sections.

Recommended concentration: For optimal performance it is recommended that the antibody be titrated for the desired application. Suggested initial dilution range:
IMC-Paraffin: 1:25 to 1:100

Description

β -actin is a 42 kDa cytoplasmic protein that is one of six different actin isoforms expressed ubiquitously in all eukaryotic cells. As an essential component of the cytoskeleton, β -actin is involved in cell motility, structure and integrity.

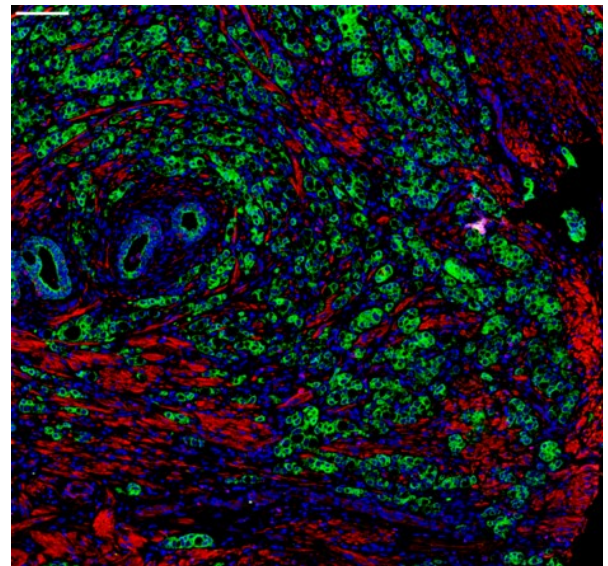

Human prostate adenocarcinoma (FFPE) stained with 154Sm-anti- β -actin (2F1-1) at a dilution of 1:50 (red pseudocolor), 148Nd-anti-pan-keratin (C51) (green pseudocolor), and iridium DNA intercalator (blue pseudocolor). Heat-mediated antigen retrieval was performed using Tris/EDTA buffer pH 9. Scale bar size = 100 μ m.

References

Chang, Q. et al. "Staining of frozen and formalin-fixed, paraffin-embedded tissues with metal-labeled antibodies for imaging mass cytometry analysis." *Current Protocols in Cytometry* 82 (2017): 12.47.1–12.47.8.

Giesen, C. et al. "Highly multiplexed imaging of tumor tissues with subcellular resolution by mass cytometry." *Nature Methods* 11 (2014): 417–22.

For technical support visit <http://techsupport.fluidigm.com>. | For general support visit www.fluidigm.com/support.

For Research Use Only. Not for use in diagnostic procedures.

This product contains antibodies manufactured by and sold under license from BioLegend[®] and licensees thereof.

Information in this publication is subject to change without notice. **Safety data sheet information:** www.fluidigm.com/sds. **Patent and license information:** www.fluidigm.com/legalnotices. **Limited Use Label License:** The purchase of this product conveys to the purchaser the limited, non-transferable right to use the purchased consumable or reagent only with Fluidigm Instruments and Systems. **EU's WEEE directive information:** www.fluidigm.com/compliance. Fluidigm, the Fluidigm logo, Hyperion, Imaging Mass Cytometry, and IMC are trademarks or registered trademarks of Fluidigm Corporation in the United States and/or other countries. All other trademarks are the sole property of their respective owners. © 2017 Fluidigm Corporation. All rights reserved. 10/2017