

Product Information Sheet

Anti-Human CD152-161Dy

Catalog #: 3161004B **Clone:** 14D3

Package Size: 100 tests Isotype: Mouse IgG2a


Storage: Store product at 4°C. Do not freeze. Formulation: Antibody stabilizer with 0.05% Sodium Azide

Cross Reactivity: Human

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF[®] analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 μ l for up to 3 X 10⁶ live cells in 100 μ l. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

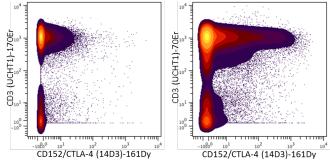

Human PBMCs were incubated for 6 hours in media alone (left) or with PMA and lonomycin (right) in the presence of monensin and brefeldin A. Cells were then fixed, permeabilized, and stained with 170Er-anti-CD3 (UCHT1) and 161Dy-anti-CD152/CTLA-4 (14D3).

Description

Cytotoxic T Lymphocyte-associated Antigen-4 (CTLA-4, CD152) is a member of the immunoglobulin superfamily and is expressed on activated T cells. CTLA-4 is a co-inhibitory molecule that functions to inhibit T cell activation by binding its ligands CD80 (B7.1) and CD86 (B7.2), which it does with higher affinity that the co-stimulatory molecule CD28. Antibodies that block the interaction of CTLA-4 with its ligands B7.1 and B7.2, including ipilimumab and tremelimumab, can enhance immune responses, including anti-tumor immunity. Recent studies have shown that CTLA-4 is expressed by a subset of regulatory T cells (Tregs).

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. Analytical Chemistry 81:6813-6822, 2009.

Ornatsky, O. I., et al. Highly multiparametric analysis by mass cytometry. J Immunol Methods 361 (1-2):1-20, 2010.

Contact Information:

Sales: sales@DVSsciences.com | Support: support@DVSsciences.com www.DVSsciences.com | For assistance by phone: 855-DVS-CYTO

This product is intended for Research Use Only. This product is not intended for therapeutic or diagnostic purposes in humans or animals.

This product contains antibodies manufactured by and sold under license from eBioscience and licensees thereof.

© 2011 DVS Sciences Inc. The trademarks mentioned herein are the property of DVS Sciences Inc. or their respective owners. DVS Sciences Inc. products are warranted to meet stated product specifications and to conform to label descriptions when used and stored properly. Unless otherwise stated, this warranty is limited to 6 months from date of sales for products used, handled and stored according to DVS Sciences, Inc. instructions. DVS Sciences, Inc. sole liability for the product is limited to replacement of the product or refund of the purchase price. Product may not be resold, modified for resale or used to manufacture commercial products without prior written approval from DVS Sciences Inc.

