

Anti-Human CD223-165Ho

Catalog: 3165028B

Package Size: 100 tests

Storage: Store product at 4°C. Do not freeze.

Reactivity: Human,

Clone: 874501


Isotype: IgG1

Formulation: Antibody stabilizer with 0.05% Sodium Azide

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF® analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 µl for up to 3×10^6 live cells in 100 µl. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

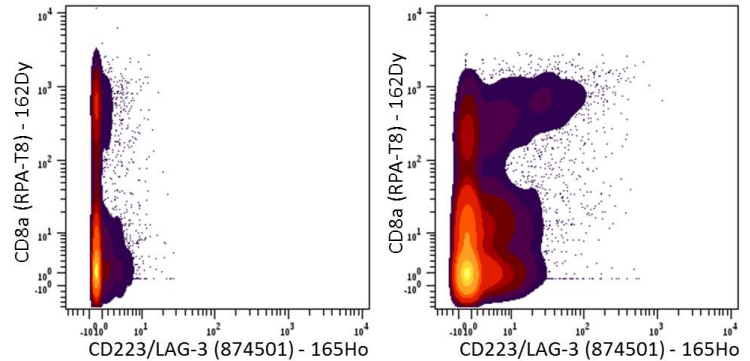

Human PBMCs were either unstimulated (left) or stimulated for 3 days with recombinant human IL-2, anti-human CD3 and anti-human CD28, and then stained with 162Dy-anti-CD8 (RPA-T8) and 165Ho-anti-CD223/LAG-3 (874501). CD3+ T cells are displayed in the analysis.

Description

Lymphocyte-activation gene 3 (LAG-3), also known as CD223, is a 70-kDa surface glycoprotein belonging to the Ig superfamily. LAG-3 shares homology with CD4, and LAG-3 binding to MHC class II is thought to contribute to the negative regulation of T cell activation. LAG-3 is expressed on T cells, including regulatory T cells and NK cells. Recently, it has been suggested that co-expression of CD49b and LAG-3 identifies human and mouse Type I regulatory T (Tr1) cells.

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. **Highly Multiparametric Analysis by Mass Cytometry.** *J Immunol Methods* 361 (1-2):1-20, 2010.

For technical support visit fluidigm.com/support

EMAIL

United States support.northamerica@fluidigm.com | Europe support.europe@fluidigm.com
 Asia techsupportasia@fluidigm.com | Latin America techsupportlatam@fluidigm.com
 All other countries techsupport@fluidigm.com

PHONE

United States (toll-free) +1 866 358 4354 | Europe +33 1 60 92 42 40
 Japan +81 3 3662 2150 | China (excluding Hong Kong) +86 21 3255 8368
 All other countries +1 650 266 6100

For Research Use Only. Not for use in diagnostic procedures.

Information in this publication is subject to change without notice. **Safety data sheet information** fluidigm.com/sds **Patent and license information** fluidigm.com/legalnotices | Fluidigm, the Fluidigm logo, Cell-ID, CyTOF, and Maxpar are trademarks or registered trademarks of Fluidigm Corporation in the United States and/or other countries. All other trademarks are the sole property of their respective owners. © 2015 Fluidigm Corporation. All rights reserved. 06/2015