

Anti-Human CD49b-161Dy

Catalog: 3161012B

Package Size: 100 tests

Storage: Store product at 4°C. Do not freeze.

Reactivity: Human,

Clone: P1E6-C5


Isotype: IgG1

Formulation: Antibody stabilizer with 0.05% Sodium Azide

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF® analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 µl for up to 3 X 10⁶ live cells in 100 µl. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

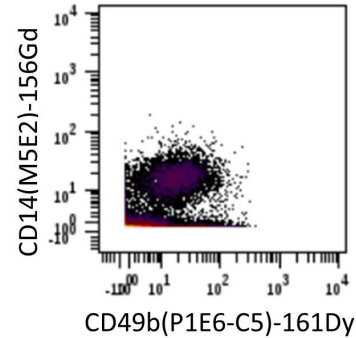

Human whole blood stained with 161Dy-anti-CD49b(P1E6-C5) and 156Gd-anti-CD14(M5E2). Total viable cells are displayed in the analysis.

Description

The P1E6-C5 monoclonal antibody specifically binds to CD49b, also known as integrin α2 chain, VLA-α2, and GP1a. CD49b is a 170kD type I transmembrane glycoprotein and belongs to the integrin family of extracellular matrix and cell-cell adhesion receptors. CD49b is expressed across a range of cell types in the blood including platelets, activated T cells, B cells, and monocytes. CD49b is reportedly associated with neonatal alloimmune thrombocytopenia. CD49b deficiencies have been associated with hemorrhagic disorders. CD49b forms the VLA-2 complex (integrin α2/β1) through a noncovalent bond with CD29. Collagen and laminin act as ligands for the VLA-2 complex. The VLA-2 complex has been reported to regulate responses mediated by proinflammatory effector cells.

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. **Highly Multiparametric Analysis by Mass Cytometry.** *J Immunol Methods* 361 (1-2):1-20, 2010.

For technical support visit fluidigm.com/support

EMAIL

United States support.northamerica@fluidigm.com | Europe support.europe@fluidigm.com
 Asia techsupportasia@fluidigm.com | Latin America techsupportlatam@fluidigm.com
 All other countries techsupport@fluidigm.com

PHONE

United States (toll-free) +1 866 358 4354 | Europe +33 1 60 92 42 40
 Japan +81 3 3662 2150 | China (excluding Hong Kong) +86 21 3255 8368
 All other countries +1 650 266 6100

For Research Use Only. Not for use in diagnostic procedures.

This product contains antibodies manufactured by and sold under license from BioLegend and licensees thereof.

Information in this publication is subject to change without notice. **Safety data sheet information** fluidigm.com/sds **Patent and license information** fluidigm.com/legalnotices | Fluidigm, the Fluidigm logo, Cell-ID, CyTOF, and Maxpar are trademarks or registered trademarks of Fluidigm Corporation in the United States and/or other countries. All other trademarks are the sole property of their respective owners. © 2015 Fluidigm Corporation. All rights reserved. 06/2015