

Anti-Human CD57-172Yb

Catalog #: 3172009B

Package Size: 100 tests

Storage: Store product at 4°C. Do not freeze.

Cross Reactivity: Human

Clone: HCD57


Isotype: Mouse Other

Formulation: Antibody stabilizer with 0.05% Sodium Azide

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF[®] analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 μ l for up to 3×10^6 live cells in 100 μ l. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

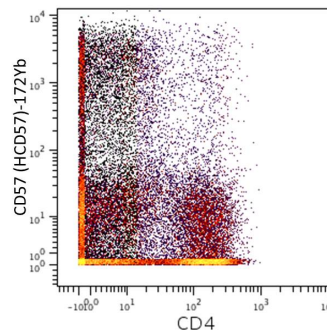

Human PBMCs stained with ¹⁴⁵Nd-anti-CD4 (RPA-T4) and ¹⁷²Yb-anti-CD57 (HCD57). Lymphocytes are displayed in the analysis.

Description

CD57, also known as HNK-1, NK-1, and Leu-7, is a 100-115 kDa cell surface glycoprotein. The function of CD57 is unknown, however it is thought to be involved in cell-matrix interactions. CD57 is expressed on a subset of peripheral blood lymphocytes including natural killer (NK) cells and CD8+ T cells. It is not expressed on red cells, granulocytes, monocytes, or platelets.

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. Highly multiparametric analysis by mass cytometry. *J Immunol Methods* 361 (1-2):1-20, 2010.

Contact Information:

Sales: sales@DVSsciences.com | Support: support@DVSsciences.com
www.DVSsciences.com | For assistance by phone: 855-DVS-CYTO

This product is intended for Research Use Only. This product is not intended for therapeutic or diagnostic purposes in humans or animals.
This product contains antibodies manufactured by and sold under license from BioLegend and licensees thereof.

© 2011 DVS Sciences Inc. The trademarks mentioned herein are the property of DVS Sciences Inc. or their respective owners. DVS Sciences Inc. products are warranted to meet stated product specifications and to conform to label descriptions when used and stored properly. Unless otherwise stated, this warranty is limited to 6 months from date of sales for products used, handled and stored according to DVS Sciences, Inc. instructions. DVS Sciences, Inc. sole liability for the product is limited to replacement of the product or refund of the purchase price. Product may not be resold, modified for resale or used to manufacture commercial products without prior written approval from DVS Sciences Inc.

