

Anti-Human HLA-DR-170Er

Catalog: 3170013B

Clone: L243

Package Size: 100 tests

Isotype: IgG2a

Storage: Store product at 4°C. Do not freeze.


Formulation: Antibody stabilizer with 0.05% Sodium Azide

Reactivity: Human

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF[®] analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 μ l for up to 3 X 10⁶ live cells in 100 μ l. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

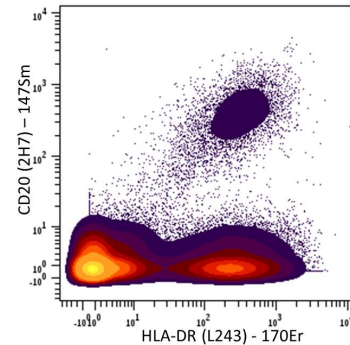

Human PBMCs stained with 147Sm-anti-CD20 (2H7) and 170Er-anti-HLA-DR (L243). Total viable cells are displayed in the analysis.

Description

HLA-DR or human leukocyte antigen DR is an MHC class II cell surface receptor that is a cell surface glycoprotein. HLA-DR is expressed on B cells, activated T cells, monocytes/macrophages, dendritic cells, and other non-professional antigen presenting cells (APCs). HLA-DR is critical for efficient antigen presentation to CD4⁺ T cells.

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Bendall, S.C., et al. Single-Cell Mass Cytometry of Differential Immune and Drug Responses Across a Human Hematopoietic Continuum. *Science* 6 May 2011: 687-696.

Bodenmiller, B. et al. State Based Profiling of Small Molecule Regulators by Cell-Multiplexed Mass Cytometry. *Nature Biotechnology* 30: 858.

Newell, E.W., et al. Cytometry by Time-of-Flight Shows Combinatorial Cytokine Expression and Virus-Specific Cell Niches within a Continuum of CD8⁺ T Cell Phenotypes. *Immunity* 36 January 2012: 142-152.

Omatsky, O. I., et al. **Highly Multiparametric Analysis by Mass Cytometry.** *J Immunol Methods* 361 (1-2):1-20, 2010.

For technical support visit fluidigm.com/support

EMAIL **United States and other countries not in Europe or Asia** techsupport@fluidigm.com **Europe** techsupport@fluidigm.com | **Asia** techsupport@fluidigm.com
PHONE **United States (toll-free)** +1 866 358 4354 | **Europe** +33 1 60 92 42 40 | **Japan** +81 3 3662 2150 | **China** +86 21 3255 8368 | **All other countries** +1 650 266 6100

For Research Use Only. Not for use in diagnostic procedures.

This product contains antibodies manufactured by and sold under license from BioLegend and licensees thereof.

Information in this publication is subject to change without notice. **Safety data sheet information** fluidigm.com/sds **Patent and license information** fluidigm.com/legalnotices | Fluidigm, the Fluidigm logo, and CyTOF are trademarks or registered trademarks of Fluidigm Corporation in the United States and/or other countries. © 2015 Fluidigm Corporation. All rights reserved. 01/2015