

Anti-Human IL-4-163Dy

Catalog: 3163011B

Package Size: 100 tests

Storage: Store product at 4°C. Do not freeze.

Reactivity: Human

Clone: MP4-25D2


Isotype: IgG1

Formulation: Antibody stabilizer with 0.05% Sodium Azide

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF[®] analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 µl for up to 3 X 10⁶ live cells in 100 µl. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

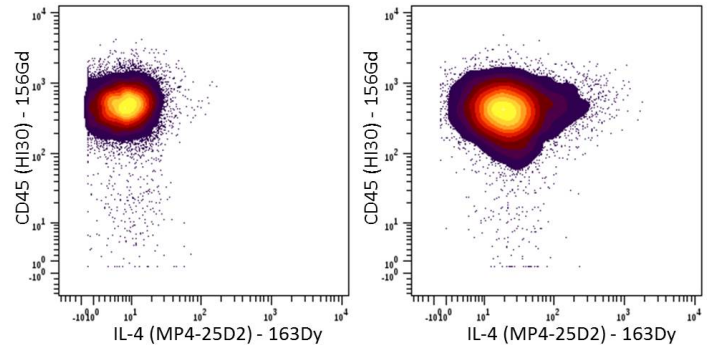

Human Th2-polarized CD4⁺ T cells were incubated for 5 hours in media alone (left) or with PMA and Ionomycin (right) in the presence of monensin and brefeldin A. Cells were then fixed, permeabilized, and stained with 156Gd-anti-CD45 (HI30) and 163Dy-anti-IL-4 (MP4-25D2).

Description

Interleukin-4 (IL-4) is a cytokine that plays a role in allergy and immunity. It is primarily produced by activated CD4⁺ T cells, mast cells and eosinophils. IL-4 stimulates TH2 cell development by activating the transcription factor STAT6, which ultimately results in enhancing expression of the cytokines IL-4, IL-5, IL-10 and IL-13. IL-4 can also regulate the production of IgE and IgG1 by B cells.

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. **Highly Multiparametric Analysis by Mass Cytometry.** *J Immunol Methods* 361 (1-2):1-20, 2010.

For technical support visit fluidigm.com/support

EMAIL **United States and other countries not in Europe or Asia** techsupport@fluidigm.com **Europe** techsupporteurope@fluidigm.com | **Asia** techsupportasia@fluidigm.com
 PHONE **United States (toll-free)** +1 866 358 4354 | **Europe** +33 1 60 92 42 40 | **Japan** +81 3 3662 2150 | **China** +86 21 3255 8368 | **All other countries** +1 650 266 6100

For Research Use Only. Not for use in diagnostic procedures.

This product contains antibodies manufactured by and sold under license from BioLegend and licensees thereof.

Information in this publication is subject to change without notice. **Safety data sheet information** fluidigm.com/sds **Patent and license information** fluidigm.com/legalnotices | Fluidigm, the Fluidigm logo, and CyTOF are trademarks or registered trademarks of Fluidigm Corporation in the United States and/or other countries. © 2015 Fluidigm Corporation. All rights reserved. 01/2015