

Anti-Human Ki-67-161Dy

Catalog #: 3161007B

Package Size: 100 tests

Storage: Store product at 4°C. Do not freeze.

Cross Reactivity: Rat, Mouse, Human, Porcine

Clone: B56


Isotype: Mouse IgG1

Formulation: Antibody stabilizer with 0.05% Sodium Azide

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF[®] analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 µl for up to 3 X 10⁶ live cells in 100 µl. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

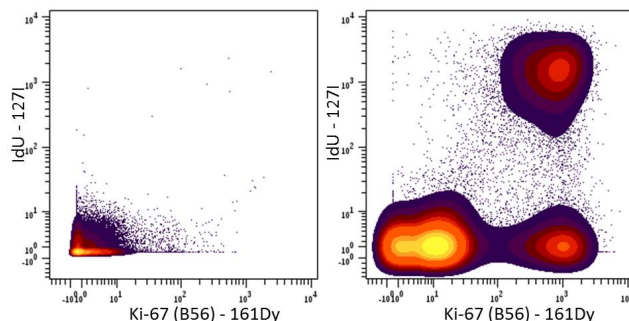

Human PBMCs were incubated for 3 days in media alone (left) or with PHA (right). Cells were then stained with IdU-1271, and then fixed, permeabilized, and stained with 161Dy-anti-Ki-67 (B56). Total viable cells are displayed in the analysis.

Description

Ki-67 protein, also known as MKI67, is a nuclear protein that is associated with cellular proliferation. Ki-67 protein is expressed in all cell types, is present during active phases of the cell cycle (G1, S, G2 and Mitosis) and not during the resting phase G0. Because of this, Ki-67 is an excellent marker for determining the fraction of proliferating cells within a given population of cells.

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Behbehani, G.K., et al. Single-cell mass cytometry adapted to measurements of the cell cycle. *Cytometry A* 81 (7): 552-566, 2012.

Ornatsky, O. I., et al. Highly multiparametric analysis by mass cytometry. *J Immunol Methods* 361 (1-2):1-20, 2010.

Contact Information:

Sales: sales@DVSsciences.com | Support: support@DVSsciences.com
www.DVSsciences.com | For assistance by phone: 855-DVS-CYTO

This product is intended for Research Use Only. This product is not intended for therapeutic or diagnostic purposes in humans or animals.

© 2011 DVS Sciences Inc. The trademarks mentioned herein are the property of DVS Sciences Inc. or their respective owners. DVS Sciences Inc. products are warranted to meet stated product specifications and to conform to label descriptions when used and stored properly. Unless otherwise stated, this warranty is limited to 6 months from date of sales for products used, handled and stored according to DVS Sciences, Inc. instructions. DVS Sciences, Inc. sole liability for the product is limited to replacement of the product or refund of the purchase price. Product may not be resold, modified for resale or used to manufacture commercial products without prior written approval from DVS Sciences Inc.

