

Anti-Human MIP-1 β -150Nd

Catalog #: 3150004B

Package Size: 100 tests

Storage: Store product at 4°C. Do not freeze.

Cross Reactivity: Human

Clone: D21-1351


Isotype: Mouse IgG1

Formulation: Antibody stabilizer with 0.05% Sodium Azide

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF[®] analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 μ l for up to 3 X 10⁶ live cells in 100 μ l. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

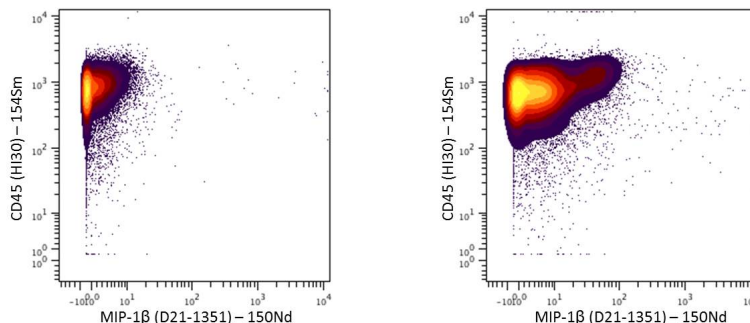

Human PBMCs were incubated for 5 hours in media alone (left) or with PMA and Ionomycin (right) in the presence of monensin and brefeldin A. Cells were then fixed, permeabilized, and stained with 154Sm-anti-CD45 (HI30) and 150Nd-anti-MIP-1 β (D21-1351).

Description

MIP-1 β (Macrophage inflammatory protein-1 β), also known as CCL4, is a cysteine-cysteine (CC) chemokine. It is a chemoattractant for natural killer cells, monocytes and a variety of other immune cells. Human MIP-1 β shares ~75% homology with mouse MIP1 β and binds to receptors, CCR5 and CCR8. The immunogen used to generate D21-1351 hybridoma was recombinant human MIP-1 β .

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Newell, E.W., et al. Cytometry by Time-of-Flight Shows Combinatorial Cytokine Expression and Virus-Specific Cell Niches within a Continuum of CD8+ T Cell Phenotypes. *Immunity* 36 January 2012: 142-152.

Ornatsky, O. I., et al. Highly multiparametric analysis by mass cytometry. *J Immunol Methods* 361 (1-2):1-20, 2010.

Contact Information:

Sales: sales@DVSciences.com | Support: support@DVSciences.com
www.DVSciences.com | For assistance by phone: 855-DVS-CYTO

This product is intended for Research Use Only. This product is not intended for therapeutic or diagnostic purposes in humans or animals.

© 2011 DVS Sciences Inc. The trademarks mentioned herein are the property of DVS Sciences Inc. or their respective owners. DVS Sciences Inc. products are warranted to meet stated product specifications and to conform to label descriptions when used and stored properly. Unless otherwise stated, this warranty is limited to 6 months from date of sales for products used, handled and stored according to DVS Sciences, Inc. instructions. DVS Sciences, Inc. sole liability for the product is limited to replacement of the product or refund of the purchase price. Product may not be resold, modified for resale or used to manufacture commercial products without prior written approval from DVS Sciences Inc.

