

Anti-Mouse Ly-6A/E (Sca-1)-169Tm

Catalog #: 3169015B

Package Size: 100 tests

Storage: Store product at 4°C. Do not freeze.

Cross Reactivity: Mouse

Clone: D7


Isotype: Rat IgG2a

Formulation: Antibody stabilizer with 0.05% Sodium Azide

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF® analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 µl for up to 3 X 10⁶ live cells in 100 µl. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

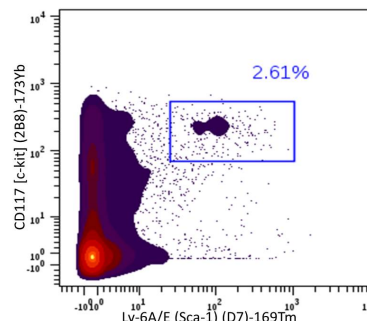

Mouse bone marrow cells stained with 173Yb-anti-CD117[c-Kit] (2B8) and 169Tm-anti-Ly-6A/E [Sca-1] (D7). Viable lineage- cells are displayed in the analysis. Lin-Sca-1+c-Kit+ LSK cells are indicated in the gate.

Description

Ly-6A/E, also known as Sca-1, is a member of the Ly-6 family of GPI-linked surface proteins. The Ly-6 family is involved in regulation and function of T cell activation. Sca-1 is a major phenotypic marker for mouse hematopoietic progenitor/stem cell subset. The D7 clone reacts with both Ly-6E.1 and Ly-6A.2 molecules expressed by mouse hematopoietic stem cells, myeloid population and peripheral T and B cells. Non-activated splenocytes from Ly-6.2 strains including C57BL, SJL, 129, AKR, and others demonstrate higher frequency of D7-positive cells compared to the Ly-6.1 strains such as BALB/c, C3H, NZB and others. Sca-1 is expressed at high levels upon activation regardless of the Ly-6 haplotype.

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. Highly multiparametric analysis by mass cytometry. *J Immunol Methods* 361 (1-2):1-20, 2010.

Contact Information:

Sales: sales@DVSsciences.com | Support: support@DVSsciences.com
www.DVSsciences.com | For assistance by phone: 855-DVS-CYTO

This product is intended for Research Use Only. This product is not intended for therapeutic or diagnostic purposes in humans or animals.
This product contains antibodies supplied by eBioscience.

© 2011 DVS Sciences Inc. The trademarks mentioned herein are the property of DVS Sciences Inc. or their respective owners. DVS Sciences Inc. products are warranted to meet stated product specifications and to conform to label descriptions when used and stored properly. Unless otherwise stated, this warranty is limited to 6 months from date of sales for products used, handled and stored according to DVS Sciences, Inc. instructions. DVS Sciences, Inc. sole liability for the product is limited to replacement of the product or refund of the purchase price. Product may not be resold, modified for resale or used to manufacture commercial products without prior written approval from DVS Sciences Inc.

